
Lentell™

OstoProtect Wipe

LENTISMED

	Pogledajte upute za uporabu / Pročitati uputstvo pre upotrebe / Preberite navodila za uporabo / Consult instructions for use		Proizvođač / Proizvođač / Proizvajalec / Manufacturer
MD	Medicinski proizvod / Medicinsko sredstvo / Medicinski pripomoček / Medical device		Datum proizvodnje / Datum proizvodnje / Datum proizvodnje / Date of manufacture
	Internetska stranica s informacijama za pacijente / Internet stranica sa informacijama za pacijente / Spletna stran z informacijami za paciente / Patient information website		Čuvati od sunčeve svjetlosti / Držati dalje od sunčeve svetlosti/ Hraniti stran od sončne svetlobe / Keep away from sunlight
LOT	Broj serije / Broj serije / Številka serije / Batch code	UDI	Jedinstveni identifikator proizvoda / Jedinstvena identifikacija proizvoda / Enotna identifikacija pripomočka / Unique device identifier
REF	Kataloški broj / Kataloški broj / Kataloška številka / Catalogue number		Distributer / Distributer / Distributer / Distributer
CE	CE oznaka / CE oznaka / Oznaka CE / CE mark		Ne koristiti ponovno / Samo za jednokratnu upotrebu / Do not re-use
	Rok uporabe / Rok upotrebe / Rok uporabe / Use-by date		Zapaljivo / Zapaljivo / Vnetljivo / Flammable

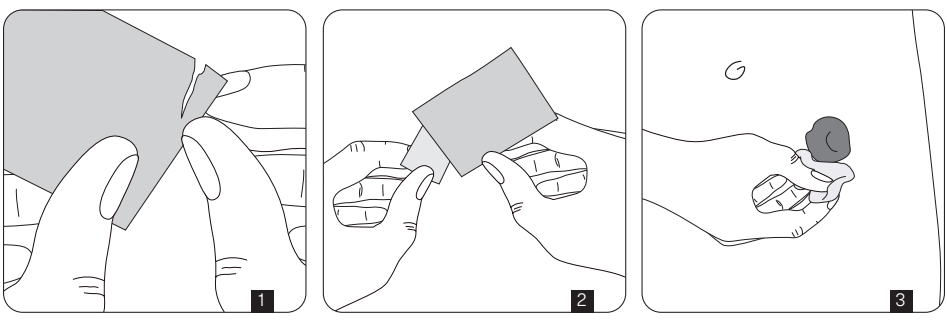HR
Upute za uporabu

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(EN)
Instructions for use

Product name, model – Lentell™ OstoProtect Wipe, barrier film wipe, ostomy accessories.
Indications for use – Lentell™ OstoProtect Wipe is intended to form a protective film on peristomal skin to help reduce irritation and support secure adhesion of ostomy appliances in patients with a colostomy, ileostomy, or urostomy, regardless of age or gender. Intended for single use by medical personnel in healthcare facilities or by patients in a home setting, on intact peristomal skin only. Please refer to the illustrated instructions for application and removal.
Storage conditions and disposal – Store the product in a cool and dry environment. Protect from sunlight. Keep away from fire and sources of ignition during storage. After use, dispose of the product according to national regulations.
Contraindications – The device must not be used in patients with known hypersensitivity to any of the product components. It must not be used in premature infants or neonates receiving care in incubators.
Precautions – This product is for external use only and must be applied to intact peristomal skin. Do not use on open wounds or damaged skin. Patch-test a small area before use. If irritation, hypersensitivity, or any adverse reaction occurs, discontinue use and consult a healthcare professional. Avoid contact with eyes and mucous membranes. In immunocompromised or otherwise clinically vulnerable patients, the device must only be used after prior consultation with a healthcare professional.
Warnings – The pouch is flammable and presents a fire hazard. Keep away from heat, sparks, and open flames during use and storage.
Safety measures – Do not use the product after the expiration date.

Any serious incident related to this product should be reported to the manufacturer and to the national competent authority.

The Lentismed logo is a registered trademark of Lentismed d.o.o. TM indicates trademarks of Lentismed d.o.o.

Lentismed d.o.o.
Slavonska Av. 26/9, 10000 Zagreb, Croatia

OPW.IFU, 1/2026
Issue date: May 2026

www.lentismed.com

T +385 1 558 4259